



QuikChange Site-Directed Mutagenesis Kit

Instruction Manual

Catalog # 200518 (30 reactions) and 200519 (10 reactions)

Revision C

Research Use Only. Not for Use in Diagnostic Procedures.

200518-12



Agilent Technologies

LIMITED PRODUCT WARRANTY

This warranty limits our liability to replacement of this product. No other warranties of any kind, express or implied, including without limitation, implied warranties of merchantability or fitness for a particular purpose, are provided by Agilent. Agilent shall have no liability for any direct, indirect, consequential, or incidental damages arising out of the use, the results of use, or the inability to use this product.

ORDERING INFORMATION AND TECHNICAL SERVICES

United States and Canada

Agilent Technologies

Stratagene Products Division

11011 North Torrey Pines Road

La Jolla, CA 92037

Telephone (858) 373-6300

Order Toll Free (800) 424-5444

Technical Services (800) 894-1304

Email techservices@agilent.com

World Wide Web www.genomics.agilent.com

Europe

Location	Telephone
Austria	01 25125 6800
Benelux	02 404 92 22
Denmark	45 70 13 00 30
Finland	010 802 220
France	0810 446 446
Germany	0800 603 1000
Italy	800 012575
Netherlands	020 547 2600
Spain	901 11 68 90
Sweden	08 506 4 8960
Switzerland	0848 8035 60
UK/Ireland	0845 712 5292

All Other Countries

Please contact your local distributor. A complete list of distributors is available at www.genomics.agilent.com.

QuikChange Site-Directed Mutagenesis Kit

CONTENTS

Materials Provided	1
Storage Conditions	1
Additional Materials Required	1
Notice to Purchaser	2
Introduction	3
QuikChange Mutagenesis Control	5
Mutagenic Primer Design	6
Primer Design Guidelines.....	6
Additional Primer Considerations	7
Protocol	8
Mutant Strand Synthesis Reaction (Thermal Cycling).....	8
<i>Dpn</i> I Digestion of the Amplification Products.....	10
Transformation of XL1-Blue Supercompetent Cells.....	10
Transformation Guidelines	12
Storage Conditions	12
Aliquoting Cells	12
Use of 14-ml BD Falcon Polypropylene Round-Bottom Tubes.....	12
Length of the Heat Pulse	12
Preparing the Agar Plates for Color Screening	12
Troubleshooting	13
Preparation of Media and Reagents	14
References	15
Endnotes	15
MSDS Information	15

QuikChange Site-Directed Mutagenesis Kit

MATERIALS PROVIDED

Materials provided	Quantity	
	Catalog #200518 ^a 30 reactions	Catalog #200519 ^b 10 reactions
<i>PfuTurbo</i> DNA polymerase (2.5 U/ μ l)	80 U	25 U
10 \times reaction buffer ^c	500 μ l	500 μ l
<i>Dpn</i> I restriction enzyme (10 U/ μ l)	300 U	100 U
Oligonucleotide control primer #1 [34-mer (100 ng/ μ l)] 5' CCA TGA TTA CGC CAA GCG CGC AAT TAA CCC TCA C 3'	750 ng	750 ng
Oligonucleotide control primer #2 [34-mer (100 ng/ μ l)] 5' GTG AGG GTT AAT TGC GCG CTT GGC GTA ATC ATG G 3'	750 ng	750 ng
pWhitescript 4.5-kb control plasmid (5 ng/ μ l)	50 ng	50 ng
dNTP mix ^{d,e}	30 μ l	10 μ l
XL1-Blue supercompetent cells ^f (blue tubes)	8 \times 200 μ l	3 \times 200 μ l
pUC18 control plasmid (0.1 ng/ μ l in TE buffer ^c)	10 μ l	10 μ l

^a The QuikChange Site-Directed Mutagenesis Kit (Catalog #200518) contains enough reagents for 30 total reactions, which includes 5 control reactions.

^b The QuikChange Site-Directed Mutagenesis Kit (Catalog #200519) contains enough reagents for 10 total reactions, which includes 5 control reactions.

^c See Preparation of Media and Reagents.

^d Thaw the dNTP mix once, prepare single-use aliquots, and store the aliquots at -20°C . **Do not subject the dNTP mix to multiple freeze-thaw cycles.**

^e **The composition of the dNTP mix is proprietary.** This reagent has been optimized for the QuikChange site-directed mutagenesis protocols and has been qualified for use in conjunction with the other kit components. Do not substitute with dNTP mixes provided with other Stratagene kits.

^f Genotype: *recA1 endA1 gyrA96 thi-1 hsdR17 supE44 relA1 lac* [F' *proAB lacI^{\Delta}ZAM15 Tn10* (Tet^r)]

STORAGE CONDITIONS

XL1-Blue Supercompetent Cells and pUC18 Control Plasmid: -80°C

All Other Components: -20°C

ADDITIONAL MATERIALS REQUIRED

14-ml BD Falcon polypropylene round-bottom tubes (BD Biosciences Catalog #352059)

5-Bromo-4-chloro-3-indolyl- β -D-galactopyranoside (X-gal)

Isopropyl- β -thio- β -D-galactopyranoside (IPTG)

NOTICE TO PURCHASER

Notice to Purchaser: Limited License

Purchase of this product includes an immunity from suit under patents specified in the product insert to use only the amount purchased for the purchaser's own internal research. No other patent rights (such as 5' Nuclease Process patent rights) are conveyed expressly, by implication, or by estoppel. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

INTRODUCTION

In vitro site-directed mutagenesis is an invaluable technique for studying protein structure-function relationships and gene expression, and for carrying out vector modification. Several approaches to this technique have been published, but these methods generally require single-stranded DNA (ssDNA) as the template¹⁻⁴ and are labor intensive or technically difficult. Our QuikChange Site-Directed Mutagenesis Kit* allows site-specific mutation in virtually any double-stranded plasmid, thus eliminating the need for subcloning into M13-based bacteriophage vectors and for ssDNA rescue.⁵ In addition, the QuikChange site-directed mutagenesis system requires no specialized vectors, unique restriction sites, or multiple transformations. This rapid four-step procedure generates mutants with greater than 80% efficiency. The protocol is simple and uses either miniprep plasmid DNA or cesium-chloride-purified DNA. For long (~8 kb) or difficult targets, we offer the QuikChange XL site directed mutagenesis kit (Catalog #200516).

The QuikChange site-directed mutagenesis kit is used to make point mutations, switch amino acids, and delete or insert single or multiple amino acids. The QuikChange site-directed mutagenesis method is performed using *PfuTurbo* DNA polymerase** and a temperature cycler. *PfuTurbo* DNA polymerase replicates both plasmid strands with high fidelity^{||} and without displacing the mutant oligonucleotide primers. The basic procedure utilizes a supercoiled double-stranded DNA (dsDNA) vector with an insert of interest and two synthetic oligonucleotide primers containing the desired mutation (see Figure 1). The oligonucleotide primers, each complementary to opposite strands of the vector, are extended during temperature cycling by *PfuTurbo* DNA polymerase. Incorporation of the oligonucleotide primers generates a mutated plasmid containing staggered nicks. Following temperature cycling, the product is treated with *Dpn* I. The *Dpn* I endonuclease (target sequence: 5'-Gm⁶ATC-3') is specific for methylated and hemimethylated DNA and is used to digest the parental DNA template and to select for mutation-containing synthesized DNA.⁶ DNA isolated from almost all *E. coli* strains is dam methylated and therefore susceptible to *Dpn* I digestion. The nicked vector DNA containing the desired mutations is then transformed into XL1-Blue supercompetent cells. The small amount of starting DNA template required to perform this method, the high fidelity of the *PfuTurbo* DNA polymerase, and the low number of thermal cycles all contribute to the high mutation efficiency and decreased potential for generating random mutations during the reaction.

Note *While plasmid DNA isolated from almost all of the commonly used E. coli strains (dam⁺) is methylated and is a suitable template for mutagenesis, plasmid DNA isolated from the exceptional dam⁻ E. coli strains, including JM110 and SCS110, is not suitable.*

* U.S. Patent Nos. 5,789,166, 5,932,419, 6,391,548, 7,132,265, 7,176,004, 5,286,632 and patents pending.

** U.S. Patent Nos. 7,045,328, 5,545,552, 5,866,395, 5,948,663, 6,183,997, 6,444,428, 6,489,150, 6,734,293 and patents pending.

^{||} *PfuTurbo* DNA polymerase has 6-fold higher fidelity in DNA synthesis than *Taq* DNA polymerase.

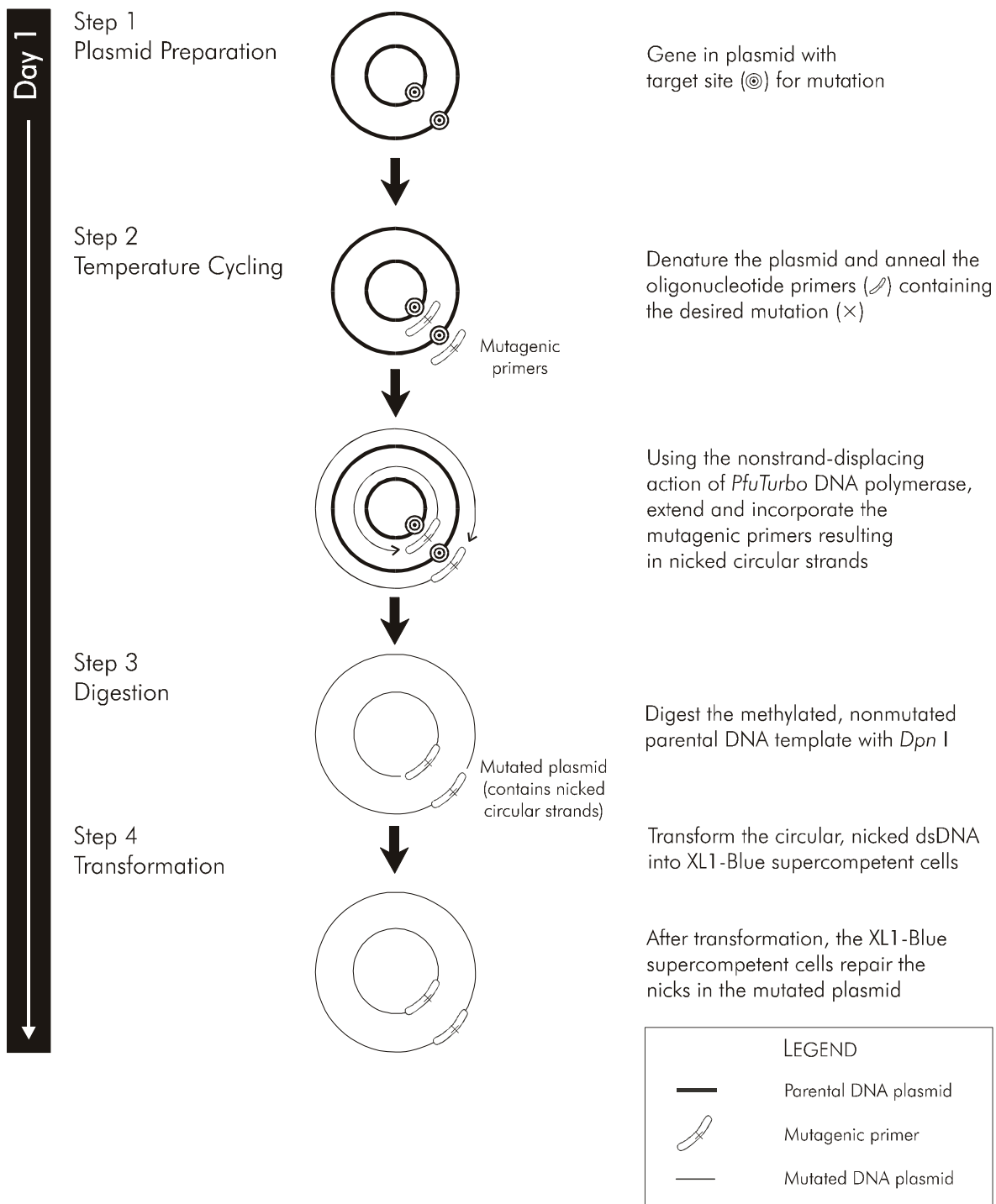


FIGURE 1 Overview of the QuikChange site-directed mutagenesis method.

QUIKCHANGE MUTAGENESIS CONTROL

The pWhitescript 4.5-kb control plasmid is used to test the efficiency of mutant plasmid generation using the QuikChange site-directed mutagenesis kit. The pWhitescript 4.5-kb control plasmid contains a stop codon (TAA) at the position where a glutamine codon (CAA) would normally appear in the β -galactosidase gene of the pBluescript II SK(–) phagemid (corresponding to amino acid 9 of the protein). XL1-Blue supercompetent cells transformed with this control plasmid appear white on LB–ampicillin agar plates (see *Preparation of Media and Reagents*), containing IPTG and X-gal, because β -galactosidase activity has been obliterated. The oligonucleotide control primers create a point mutation on the pWhitescript 4.5-kb control plasmid that reverts the T residue of the stop codon (TAA) at amino acid 9 of the β -galactosidase gene to a C residue, to produce the glutamine codon (CAA) found in the wild-type sequence. Following transformation, colonies can be screened for the β -galactosidase (β -gal⁺, blue) phenotype.

MUTAGENIC PRIMER DESIGN

Note *Mutagenic primers can be designed using our web-based QuikChange Primer Design Program available online at www.agilent.com/genomics/qcpd.*

Primer Design Guidelines

The mutagenic oligonucleotide primers for use in this protocol must be designed individually according to the desired mutation. The following considerations should be made for designing mutagenic primers:

- ◆ Both of the mutagenic primers must contain the desired mutation and anneal to the same sequence on opposite strands of the plasmid.
- ◆ Primers should be between 25 and 45 bases in length, with a melting temperature (T_m) of $\geq 78^\circ\text{C}$. Primers longer than 45 bases may be used, but using longer primers increases the likelihood of secondary structure formation, which may affect the efficiency of the mutagenesis reaction. The following formula is commonly used for estimating the T_m of primers:

$$T_m = 81.5 + 0.41(\%GC) - (675/N) - \% \text{ mismatch}$$

For calculating T_m :

- N is the primer length in bases
- values for **%GC** and **% mismatch** are whole numbers

For calculating T_m for primers intended to introduce insertions or deletions, use this modified version of the above formula:

$$T_m = 81.5 + 0.41(\%GC) - (675/N)$$

where N does not include the bases which are being inserted or deleted.

Note *When using primer design software for QuikChange site-directed mutagenesis applications, be aware that the T_m calculated by the primer design software may differ from the T_m value calculated using the formula presented above. We recommend verifying primer T_m 's using the formula above or by using the QuikChange T_m calculator, available online at www.genomics.agilent.com.*

- ◆ The desired mutation (deletion or insertion) should be in the middle of the primer with ~10–15 bases of correct sequence on both sides.
- ◆ The primers optimally should have a minimum GC content of 40% and should terminate in one or more C or G bases.

Additional Primer Considerations

- ◆ The mutagenesis protocol uses 125 ng of each oligonucleotide primer. To convert nanograms to picomoles of oligo, use the following equation:

$$X \text{ pmoles of oligo} = \frac{\text{ng of oligo}}{330 \times \# \text{ of bases in oligo}} \times 1000$$

For example, for 125 ng of a 25-mer:

$$\frac{125 \text{ ng of oligo}}{330 \times 25 \text{ bases}} \times 1000 = 15 \text{ pmole}$$

- ◆ Primers need not be 5' phosphorylated but **must** be purified either by fast polynucleotide liquid chromatography (FPLC) or by polyacrylamide gel electrophoresis (PAGE). Failure to purify the primers results in a significant decrease in mutation efficiency.
- ◆ It is important to keep primer concentration in excess. We suggest varying the amount of template while keeping the concentration of the primer constantly in excess.

PROTOCOL

Mutant Strand Synthesis Reaction (Thermal Cycling)

Notes *Ensure that the plasmid DNA template is isolated from a dam^+ E. coli strain. The majority of the commonly used E. coli strains are dam^+ . Plasmid DNA isolated from dam^- strains (e.g. JM110 and SCS110) is not suitable.*

To maximize temperature-cycling performance, We strongly recommend using thin-walled tubes, which ensure ideal contact with the temperature cycler's heat blocks. The following protocols were optimized using thin-walled tubes.

1. Synthesize two complimentary oligonucleotides containing the desired mutation, flanked by unmodified nucleotide sequence. Purify these oligonucleotide "primers" prior to use in the following steps (see *Mutagenic Primer Design*).

2. Prepare the control reaction as indicated below:

5 μl of 10 \times reaction buffer (see *Preparation of Media and Reagents*)
2 μl (10 ng) of pWhitescript 4.5-kb control plasmid (5 ng/ μl)
1.25 μl (125 ng) of oligonucleotide control primer #1
[34-mer (100 ng/ μl)]
1.25 μl (125 ng) of oligonucleotide control primer #2
[34-mer (100 ng/ μl)]
1 μl of dNTP mix
39.5 μl of double-distilled water (ddH_2O) to a final volume of 50 μl

Then add

1 μl of *PfuTurbo* DNA polymerase (2.5 U/ μl)

3. Prepare the sample reaction(s) as indicated below:

Note *Set up a series of sample reactions using various concentrations of dsDNA template ranging from 5 to 50 ng (e.g., 5, 10, 20, and 50 ng of dsDNA template) while keeping the primer concentration constant.*

5 μl of 10 \times reaction buffer
 X μl (5–50 ng) of dsDNA template
 X μl (125 ng) of oligonucleotide primer #1
 X μl (125 ng) of oligonucleotide primer #2
1 μl of dNTP mix
 ddH_2O to a final volume of 50 μl

Then add

1 μl of *PfuTurbo* DNA polymerase (2.5 U/ μl)

- If the thermal cycler to be used does not have a hot-top assembly, overlay each reaction with ~30 μ l of mineral oil.

TABLE I

Cycling Parameters for the QuikChange Site-Directed Mutagenesis Method

Segment	Cycles	Temperature	Time
1	1	95°C	30 seconds
2	12–18	95°C	30 seconds
		55°C	1 minute
		68°C	1 minute/kb of plasmid length*

* For example, a 5-kb plasmid requires 5 minutes at 68°C per cycle.

- Cycle each reaction using the cycling parameters outlined in Table I. (For the control reaction, use a 5-minute extension time and run the reaction for 18 cycles.)
- Adjust segment 2 of the cycling parameters in accordance with the type of mutation desired (see the following table):

Type of mutation desired	Number of cycles
Point mutations	12
Single amino acid changes	16
Multiple amino acid deletions or insertions	18

- Following temperature cycling, place the reaction on ice for 2 minutes to cool the reaction to $\leq 37^\circ\text{C}$.

Note *If desired, amplification may be checked by electrophoresis of 10 μ l of the product on a 1% agarose gel. A band may or may not be visualized at this stage. In either case proceed with Dpn I digestion and transformation.*

Dpn I Digestion of the Amplification Products

Note *It is important to insert the pipet tip below the mineral oil overlay when adding the Dpn I restriction enzyme to the reaction tubes during the digestion step or when transferring the 1 µl of Dpn I-treated DNA to the transformation reaction. Using specialized aerosol-resistant pipet tips, which are small and pointed, facilitates this process.*

1. Add 1 µl of the Dpn I restriction enzyme (10 U/µl) directly to each amplification reaction below the mineral oil overlay using a small, pointed pipet tip.
2. Gently and thoroughly mix each reaction mixture by pipetting the solution up and down several times. Spin down the reaction mixtures in a microcentrifuge for 1 minute and immediately incubate each reaction at 37°C for 1 hour to digest the parental (i.e., the nonmutated) supercoiled dsDNA.

Transformation of XL1-Blue Supercompetent Cells

Notes *Please read the Transformation Guidelines before proceeding with the transformation protocol.*

XL1-Blue cells are resistant to tetracycline. If the mutagenized plasmid contains only the tet^R resistance marker, an alternative tetracycline-sensitive strain of competent cells must be used.

1. Gently thaw the XL1-Blue supercompetent cells on ice. For each control and sample reaction to be transformed, aliquot 50 µl of the supercompetent cells to a *prechilled* 14-ml BD Falcon polypropylene round-bottom tube.
2. Transfer 1 µl of the Dpn I-treated DNA from each control and sample reaction to separate aliquots of the supercompetent cells.

Note *Carefully remove any residual mineral oil from the pipet tip before transferring the Dpn I-treated DNA to the transformation reaction.*

As an optional control, verify the transformation efficiency of the XL1-Blue supercompetent cells by adding 1 µl of the pUC18 control plasmid (0.1 ng/µl) to a 50-µl aliquot of the supercompetent cells.

Swirl the transformation reactions gently to mix and incubate the reactions on ice for 30 minutes.

3. Heat pulse the transformation reactions for 45 seconds at 42°C and then place the reactions on ice for 2 minutes.

Note *This heat pulse has been optimized for transformation in 14-ml BD Falcon polypropylene round-bottom tubes.*

4. Add 0.5 ml of NZY⁺ broth (see *Preparation of Media and Reagents*) preheated to 42°C and incubate the transformation reactions at 37°C for 1 hour with shaking at 225–250 rpm.
5. Plate the appropriate volume of each transformation reaction, as indicated in the table below, on agar plates containing the appropriate antibiotic for the plasmid vector.

For the mutagenesis and transformation controls, spread cells on LB–ampicillin agar plates containing 80 µg/ml X-gal and 20 mM IPTG (see *Preparing the Agar Plates for Color Screening*).

Transformation reaction plating volumes

Reaction Type	Volume to Plate
pWhitescript mutagenesis control	250 µl
pUC18 transformation control	5 µl (in 200 µl of NZY ⁺ broth)*
Sample mutagenesis	250 µl on each of two plates (entire transformation reaction)

* Place a 200-µl pool of NZY⁺ broth on the agar plate, pipet the 5 µl of the transformation reaction into the pool, then spread the mixture.

6. Incubate the transformation plates at 37°C for >16 hours.

Expected Results for the Control Transformations

The expected colony number from the transformation of the pWhitescript control mutagenesis reaction is between 50 and 800 colonies. Greater than 80% of the colonies should contain the mutation and appear as blue colonies on agar plates containing IPTG and X-gal.

Note *The mutagenesis efficiency (ME) for the pWhitescript 4.5-kb control plasmid is calculated by the following formula:*

$$ME = \frac{\text{Number of blue colony forming units (cfu)}}{\text{Total number of colony forming units (cfu)}} \times 100\%$$

If transformation of the pUC18 control plasmid was performed, >250 colonies should be observed (transformation efficiency >10⁸ cfu/µg) with >98% of the colonies having the blue phenotype.

Expected Results for Sample Transformations

The expected colony number is between 10 and 1000 colonies, depending upon the base composition and length of the DNA template employed. For suggestions on increasing colony number, see *Troubleshooting*. The insert of interest should be sequenced to verify that selected clones contain the desired mutation(s).

TRANSFORMATION GUIDELINES

It is important to store the XL1-Blue supercompetent cells at -80°C to prevent a loss of efficiency. For best results, please follow the directions outlined in the following sections.

Storage Conditions

The XL1-Blue supercompetent cells are very sensitive to even small variations in temperature and must be stored at the bottom of a -80°C freezer. Transferring tubes from one freezer to another may result in a loss of efficiency. The XL1-Blue supercompetent cells should be placed at -80°C directly from the dry ice shipping container.

Aliquoting Cells

When aliquoting, keep the XL1-Blue supercompetent cells on ice at all times. It is essential that the BD Falcon polypropylene tubes are placed on ice before the cells are thawed and that the cells are aliquoted directly into the prechilled tubes.

Use of 14-ml BD Falcon Polypropylene Round-Bottom Tubes

It is important that 14-ml BD Falcon polypropylene round-bottom tubes (BD Biosciences Catalog #352059) are used for the transformation protocol because the duration of the heat-pulse step is critical and has been optimized for the thickness and shape of these tubes.

Length of the Heat Pulse

There is a defined "window" of highest efficiency for the XL1-Blue supercompetent cells resulting from the heat pulse in step 3 of the transformation protocol. Optimal efficiencies are observed when cells are heat pulsed for 45 seconds. Heat pulsing for at least 45 seconds is recommended to allow for slight variations in the length of incubation. Efficiencies decrease sharply when pulsing for <30 seconds or for >45 seconds.

Preparing the Agar Plates for Color Screening

To prepare the LB agar plates for blue-white color screening, add 80 $\mu\text{g}/\text{ml}$ of 5-bromo-4-chloro-3-indolyl- β -D-galactopyranoside (X-gal), 20 mM isopropyl-1-thio- β -D-galactopyranoside (IPTG), and the appropriate antibiotic to the LB agar. Alternatively, 100 μl of 10 mM IPTG and 100 μl of 2% X-gal can be spread on the LB agar plates 30 minutes prior to plating the transformations. Prepare the IPTG in sterile dH_2O ; prepare the X-gal in dimethylformamide (DMF). Do not mix the IPTG and the X-gal before pipetting them onto the plates because these chemicals may precipitate.

TROUBLESHOOTING

When used according to the guidelines outlined in this instruction manual, this kit provides a reliable means to conduct site-directed mutagenesis using dsDNA templates. Variations in the base composition and length of the DNA template and in thermal cycler performance may contribute to differences in mutagenesis efficiency. We provide the following guidelines for troubleshooting these variations.

Observation	Suggestion(s)
Low transformation efficiency or low colony number	Ensure that excess mineral oil is not transferred into the transformation reaction when pipetting the <i>Dpn</i> I-treated DNA. Using the smallest pipet tips available, insert the pipet tip completely below the mineral layer overlay and clear the pipet tip while submerged beneath the mineral oil overlay before collecting the sample.
	Ensure that sufficient mutant DNA is synthesized in the reaction. Increase the amount of the <i>Dpn</i> I-treated DNA used in the transformation reaction to 4 μ l.
	Visualize the DNA template on a gel to verify the quantity and quality. Nicked or linearized plasmid DNA will not generate complete circular product. Verify that the template DNA is at least 80% supercoiled.
	It is not uncommon to observe low numbers of colonies, especially when generating large mutations. Most of the colonies that do appear, however, will contain mutagenized plasmid.
	Ethanol precipitate the <i>Dpn</i> I digested PCR product, and resuspend in a decreased volume of water before transformation.
Low mutagenesis efficiency or low colony number with the control reaction	Different thermal cyclers may contribute to variations in ramping efficiencies. Adjust the cycling parameters for the control reaction and repeat the protocol for the sample reactions.
	Ensure that supercompetent cells are stored at the bottom of a -80°C freezer immediately upon arrival (see also <i>Transformation Guidelines</i>).
	Verify that the agar plates were prepared correctly. See <i>Preparing the Agar Plates for Color Screening</i> , and follow the recommendations for IPTG and X-Gal concentrations carefully.
	For best visualization of the blue (β -gal ⁺) phenotype, the control plates must be incubated for at least 16 hours at 37°C .
	Avoid multiple freeze-thaw cycles for the dNTP mix. Thaw the dNTP mix once, prepare single-use aliquots, and store the aliquots at -20°C . Do not subject the dNTP mix to multiple freeze-thaw cycles.
	Adjust the cycling parameters for the sample reaction to overcome differences in ramping efficiencies of thermal cyclers.
Low mutagenesis efficiency with the sample reaction(s)	Add the <i>Dpn</i> I restriction enzyme below the mineral oil overlay in the digestion step and ensure proper mixing of all components in the reaction especially the <i>Dpn</i> I.
	Allow sufficient time for the <i>Dpn</i> I to completely digest the parental template; repeat the digestion if too much DNA template was present.
	Avoid multiple freeze-thaw cycles for the dNTP mix. Thaw the dNTP mix once, prepare single-use aliquots, and store the aliquots at -20°C . Do not subject the dNTP mix to multiple freeze-thaw cycles.
	The formation of secondary structures may be inhibiting the mutagenesis reaction. Increasing the annealing temperature up to 68°C may help to alleviate secondary structure formation and improve mutagenesis efficiency.

Table continues on the following page

Table continues from the previous page

False positives	Poor quality primers can lead to false positives. Radiolabel the primers and check for degradation on an acrylamide gel or resynthesize the primers.
	False priming can lead to false positives. Increase the stringency of the reaction by increasing the annealing temperature up to 68°C.
Unwanted deletion or recombination of plasmid DNA following mutagenesis and transformation	Transform the mutagenesis reaction into competent cells that are designed to prevent recombination events, such as Stratagene SURE 2 Supercompetent Cells (Catalog #200152). Note that SURE 2 competent cells are not recommended for use with mutagenized plasmids greater than 10 kb in size; note also that SURE 2 cells are Kan ^r , Tet ^r , and Chl ^r , and are not compatible with plasmid selection using kanamycin, tetracycline, or chloramphenicol resistance markers.

PREPARATION OF MEDIA AND REAGENTS

LB Agar (per Liter) 10 g of NaCl 10 g of tryptone 5 g of yeast extract 20 g of agar Add deionized H ₂ O to a final volume of 1 liter Adjust pH to 7.0 with 5 N NaOH Autoclave Pour into petri dishes (~25 ml/100-mm plate)	LB–Ampicillin Agar (per Liter) 1 liter of LB agar, autoclaved Cool to 55°C Add 10 ml of 10-mg/ml filter-sterilized ampicillin Pour into petri dishes (~25 ml/100-mm plate)
NZY⁺ Broth (per Liter) 10 g of NZ amine (casein hydrolysate) 5 g of yeast extract 5 g of NaCl Add deionized H ₂ O to a final volume of 1 liter Adjust to pH 7.5 using NaOH Autoclave Add the following filter-sterilized supplements prior to use: 12.5 ml of 1 M MgCl ₂ 12.5 ml of 1 M MgSO ₄ 20 ml of 20% (w/v) glucose (or 10 ml of 2 M glucose)	10× Reaction Buffer 100 mM KCl 100 mM (NH ₄) ₂ SO ₄ 200 mM Tris-HCl (pH 8.8) 20 mM MgSO ₄ 1% Triton [®] X-100 1 mg/ml nuclease-free bovine serum albumin (BSA)
	TE Buffer 10 mM Tris-HCl (pH 7.5) 1 mM EDTA

REFERENCES

1. Kunkel, T. A. (1985) *Proc Natl Acad Sci U S A* 82(2):488-92.
2. Sugimoto, M., Esaki, N., Tanaka, H. and Soda, K. (1989) *Anal Biochem* 179(2):309-11.
3. Taylor, J. W., Ott, J. and Eckstein, F. (1985) *Nucleic Acids Res* 13(24):8765-85.
4. Vandeyar, M. A., Weiner, M. P., Hutton, C. J. and Batt, C. A. (1988) *Gene* 65(1):129-33.
5. Papworth, C., Bauer, J. C., Braman, J. and Wright, D. A. (1996) *Strategies* 9(3):3-4.
6. Nelson, M. and McClelland, M. (1992) *Methods Enzymol* 216:279-303.

ENDNOTES

Triton® is a registered trademark of Rohm and Haas Co.

MSDS INFORMATION

The Material Safety Data Sheet (MSDS) information for Stratagene products is provided on the web at <http://www.genomics.agilent.com>. MSDS documents are not included with product shipments.

STRATAGENE

An Agilent Technologies Division

QuikChange Site-Directed Mutagenesis Kit

Catalog #200518 and #200519

QUICK-REFERENCE PROTOCOL

- ♦ Prepare the control and sample reaction(s) as indicated below:

Note Set up a series of sample reactions using various concentrations ranging from 5 to 50 ng of dsDNA template (e.g., 5, 10, 20, and 50 ng of dsDNA template).

Control Reaction

5 μ l of 10 \times reaction buffer
2 μ l (10 ng) of pWhitescript 4.5-kb control template (5 ng/ μ l)
1.25 μ l (125 ng) of oligonucleotide control primer #1 [34-mer (100 ng/ μ l)]
1.25 μ l (125 ng) of oligonucleotide control primer #2 [34-mer (100 ng/ μ l)]
1 μ l of dNTP mix
ddH₂O to a final volume of 50 μ l

Sample Reaction

5 μ l of 10 \times reaction buffer
X μ l (5–50 ng) of dsDNA template
X μ l (125 ng) of oligonucleotide primer #1
X μ l (125 ng) of oligonucleotide primer #2
1 μ l of dNTP mix
ddH₂O to a final volume of 50 μ l

- ♦ Then add 1 μ l of *PfuTurbo* DNA polymerase (2.5 U/ μ l) to each control and sample reaction.
- ♦ Overlay each reaction with 30 μ l of mineral oil.
- ♦ Cycle each reaction using the cycling parameters outlined in the following table:

Segment	Cycles	Temperature	Time
1	1	95°C	30 seconds
2	12–18	95°C	30 seconds
		55°C	1 minute
		68°C	1 minute/kb of plasmid length

- ♦ Adjust segment 2 of the cycling parameters in accordance with the type of mutation desired (see the table in step 6 of *Mutant Strand Synthesis Reaction (Thermal Cycling)* in the instruction manual).
- ♦ Add 1 μ l of the *Dpn* I restriction enzyme (10 U/ μ l) below the mineral oil overlay.
- ♦ Gently and thoroughly mix each reaction, spin down in a microcentrifuge for 1 minute, and immediately incubate at 37°C for 1 hour to digest the parental supercoiled dsDNA.
- ♦ Transform 1 μ l of the *Dpn* I-treated DNA from each control and sample reaction into separate 50- μ l aliquots of XL1-Blue supercompetent cells (see *Transformation of XL1-Blue Supercompetent Cells* in the instruction manual).